

Simple steps for a better experience



Get the best from your Grenco Medical Elite II
— and keep it performing at its best.

1

Medium grind, packed loose —
don't overfill the chamber.

2

Slow, steady inhales work best.

3

Start low, work up to your preferred temperature.
Sessions typically last 3–5 min.

4

Let cool, then use the built-in pick tool to
empty & repack.

Keeping your device in top condition.

PLEASE AVOID

- Submerging your device in liquids.
- Abrasive cleaners
(*only use Isopropyl Alcohol*)
or cleaning while powered on
- Leaving within reach of children



Scan for the
in-depth guide



SPONSOR AUSTRALIA

MS MOOD STATE
MEDICAL

MOOD STATE MEDICAL
moodstatemedical.au | 1300 069 385
support@moodstatemedical.au
Blaxland House Suite 6/5-7 Ross Street
Parramatta NSW 2150 AUSTRALIA

ARTG# 526764
Class IIb medical device

Quick Start Guide



Register

- Register your device at GrencoMedical.com/Register




Charge

- Use included USB-C cable
- Charge fully before first use



Load

- Remove mouthpiece
- Fill chamber with ground material
- Re-attach mouthpiece
- Pack with pick tool (do not overpack)

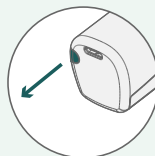

Pick Tool

- Stored at the base of device — use for packing, stirring & clearing the chamber



Activate

Hold for 3 seconds to power on / off

+ / - to set desired temperature

Tap 2X to begin / cancel a heat session

Green LED on - Temperature has been reached

Tap 5X to access device settings menu



Wait

- Wait 30 seconds after device reaches temperature & vibrates for the flower to heat before inhaling



Inhale

- LED on — inhale from mouthpiece - Temperature has been reached

Clean

- Clean screen & chamber after every use
- Deep clean: remove mouthpiece insert, wipe with isopropyl alcohol

**Dry all parts before reassembling*

Settings

Tap 5X to access device settings menu
Navigate with + / - , select with power

- Screen brightness
- Haptic feedback on / off
- Temperature unit °F / °C