




grenco
medical



GRENCO MEDICAL ELITE II
INSTRUCTIONS FOR USE

GRENCO MEDICAL ELITE II EXPLODED VIEW

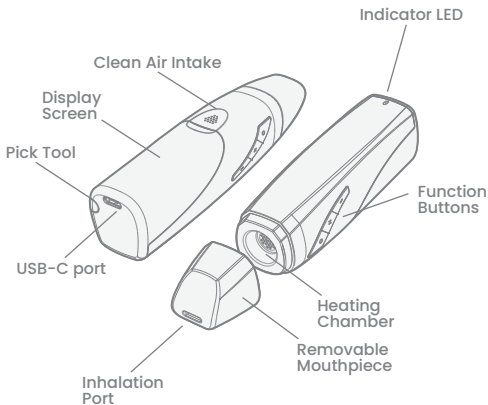

Throughout the expected lifetime of the medical device, certain components may require replacement. Should any part of the medical device require replacement, please contact our customer support team. For assistance, contact GS HOLISTIC, LLC at SUPPORT@GRENCOMEDICAL.COM or call +1 833-691-3224.



EXPLANATION OF PACKAGING AND INSTRUCTION MANUAL SYMBOLS



Read instructions symbol. Please read and follow these instructions carefully.



Warning symbol. To prevent harm to individuals and protect medical device, compliance with instructions featuring this symbol is compulsory.



Safety symbol. To prevent medical device damage, it is essential to adhere to the instructions indicated by this symbol.



Manufacturer symbol. The manufacturer's name and address are located adjacent to this symbol.

SN

Serial number symbol. Following this symbol is the serial number of the medical device.

LOT

Lot number symbol. Following this symbol is the lot number of the medical device production.



Manufacturing date symbol. Following this symbol is the date that the medical device was manufactured.



Expiry date symbol. This symbol will be next to the date (year/ month) by which the medical device must be used by.



Do not use the medical device if the box is damaged.



Protect medical device from moisture and humidity.



Protect medical device from direct sunlight.



These are Applied Parts that make direct electrical contact with the patient.



The medical device should not be discarded as unsorted waste but must be sent to separate collection facilities for recovery and recycling.

REF

Device Name: GRENCO MEDICAL ELITE II





WARNINGS

- The medical device is solely intended for vaporising the recommended cannabinoids found in dried flowers.
- Individuals requiring assistance must be under constant supervision during application as the medical device contains small components that could obstruct the airway and potentially cause choking. Therefore, always ensure that the medical device and its accessories are kept out of reach of infants and children.
- When charging the medical device battery, exclusively use the provided USB-C cable for the medical device and avoid unsupervised charging.
- The medical device offers a highly efficient and safe method for alveolar cannabinoid application.
- Before commencing treatment, patients must undergo a medical evaluation by a licensed physician.
- Usage is reserved for patients aged 18 or older, after thoroughly reviewing these instructions and following the physician's recommendations.
- Adults can operate the medical device by adhering to these Instructions for Use or as directed by their physician.

- Refrain from using the medical device if the patient has respiratory or lung conditions, as vapor exposure may lead to irritation and coughing.
- Only employ medicinal cannabis flowers in the medical device as per your physician's prescription.
- The ideal temperature range for cannabinoid extraction falls between 180°C and 210°C. Once inhaled, cannabinoids are rapidly absorbed in the alveoli, and their effects can be experienced within minutes.

STORAGE AND TRANSPORT CONDITIONS





SAFETY AND CAUTIONS

- Inspect the packaging and housing of the medical device. If you notice any damage to the medical device housing, refrain from using the medical device and promptly get in touch with customer support.
- When recharging the battery, make sure to exclusively utilize the provided USB-C cable that comes with the medical device.
- Clean the medical device solely with a dry cloth, or if necessary, a damp cloth. Never submerge the medical device or the power adapter in any liquids.
- Avoid exposing the medical device to direct sunlight, rain, or extreme temperatures. Drastic fluctuations in ambient temperature can lead to condensation and impair the medical device's functionality.
- Operate the medical device within an ambient temperature range of +5°C to +30°C, maintaining an ambient air relative humidity between 15% and 90%, and ensuring an atmospheric pressure between 700 hPa and 1050 hPa. If the mouthpiece temperature exceeds +43°C, the medical device will auto shut off for 20 minutes to allow mouthpiece to cool to safe operating temperature.

INTENDED USE

THE GRENCO MEDICAL ELITE II VAPORISER IS INTENDED FOR USE AT HOME, IN A HOSPITAL, OR DOCTOR'S OFFICE FOR VAPORISING OF CANNABIS FLOWERS AND FOR INHALING CANNABINOIDS DERIVED THEREFROM. ACTIVE INGREDIENTS IN THE INHALED CANNABINOIDS ARE DELIVERED INTO THE HUMAN BODY VIA THE ALVEOLI FOR THE APPROPRIATE MEDICAL INDICATION AS PRESCRIBED BY YOUR PHYSICIAN. USE IS INTENDED BY PATIENTS WHO ARE OLDER THAN 18 YEARS OF AGE.

CONTRAINDICATIONS

- The medical device does not have any device-specific contraindications.
- Nevertheless, if the user has a pre-existing respiratory tract or lung condition, the medical device should only be employed following consultation with and approval from their attending physician. The vapors may induce irritation in the respiratory tract and lungs, potentially resulting in bouts of coughing, depending on vapor density. can be experienced within minutes.



UNPACKING AND COMMISSIONING OF THE MEDICAL DEVICE

- Remove the medical device from the box and all the accessories.
- Register serial number found at the bottom of medical device or packaging at grencomedical.com/register.

- Fully charge medical device using included USB-C Charging Cable.

- Grind medicinal cannabis flowers.

- Remove mouthpiece and load the ground medicinal cannabis flowers into the medical device heating chamber.

- Re-attach mouthpiece.

- Press and hold power button for three seconds to power on/off the medical device.

- Use + and - buttons to set desired temperature.

- Double-click power button to begin/cancel heat session.

- When indicator LED turns green, inhale from the mouthpiece.



DEVICE SETTINGS

- Press power button five times to access medical device Settings Menu.
- Press + and - buttons to navigate menu. Press power button to select: Screen brightness, Toggling haptic Feedback, Set displaying between °F/°C.
- Access pick tool at bottom of device for packing ground material into heating chamber, stirring of ground material, and removing of spent material from the heating chamber.



CONTACT AND SUPPORT

GS HOLISTIC, LLC
1123 LOCUST ST. #166
ST. LOUIS, MO 63101 USA

SUPPORT@GRENCOMEDICAL.COM
+1 833-691-3224



MANUFACTURED BY
DONGGUAN MEDICA TECHNOLOGY CO., LTD
ROOM 408, BUILDING NO. 1
NO. 136 GANG JIAN ROAD
CHANG PING TOWN, DONGGUAN
GUANGDONG 523581 CHINA

SPONSOR AUSTRALIA



MOOD STATE MEDICAL
moodstatemedical.au | 1300 069 385
support@moodstatemedical.au
Blaxland House Suite 6/5-7 Ross Street
Parramatta NSW 2150 AUSTRALIA